

HAIYE



HYHJ-KC Reusable Video Laryngoscope



Application


It is provided for clinical use to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and implement mechanical ventilation.

Suitable for anesthesia departments, ICU,



Feature and Configuration

The host display and camera parts are connected by plug and pull mode, which is efficient, strong and beautiful. The product can be used with 6 different specifications of blades can be reused after disinfection, its specifications are Miller 0, Miller 1, Mac 2, Mac 3, Mac 4, D-blade, can be applied to all people.

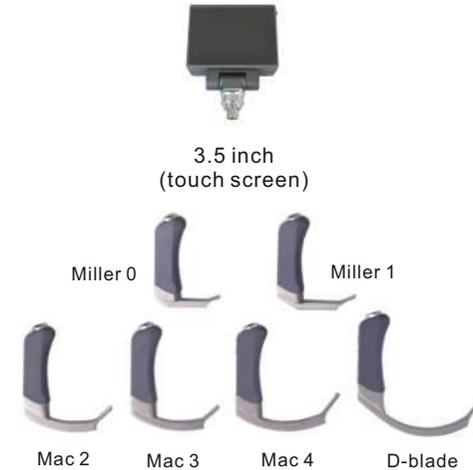

Parameter

Monitor	3" LCD display / 3.5" color touch screen display
Camera	pixel: ≥3million; light illuminance: ≥3000 lux
Battery	battery life: ≥4h; charging time: ≤5h; capacity: 3400mAh
Interface Function	HDMI, Type-C; Photo, video, freeze, playback, store and export data



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HYHJ-1320 Disposable Video Laryngoscope



Application

It is provided for clinical use to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and implement mechanical ventilation.

Suitable for anesthesia departments, ICU, pediatrics, emergency centers, ambulances, etc. in medical institutions.



Feature and Configuration

The product adopts the structural design combining stainless steel and polymer material, with light weight, high strength, full functions, long life; the product is equipped with camera of 2 million pixels, 3" color display screen, standard with three specifications of camera system, corresponding to three specifications of disposable laryngoscopes blades, can be applied to all people.



Camera system



Single-use laryngoscope blade

Parameter

Monitor	3" color display
Camera	Pixel: ≥ 2 million; illuminance: ≥ 3000 lux
Battery	Battery life: ≥ 6 h; charging time: ≤ 4 h; capacity: 3400mAh
Function	Take pictures, video, freeze, store and export data
Storage	-40°C~70°C ambient temperature
Humidity	$\leq 93\%$ rh



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com



HYHJ-1330 Disposable Video laryngoscope

/HYHJ-1330 I, HYHJ-1330 II, HYHJ-1330 III/



Application

It is provided for clinical use to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and implement mechanical ventilation.

Suitable for anesthesia departments, ICU, pediatrics, emergency centers, ambulances, etc. in medical institutions.



Feature and Configuration

This product has a simple and durable structure, and is designed as an all-in-one machine with one specification, applicable to Adults, children, newborns; corresponding to 3 kinds of disposable laryngoscope blade, can be applied to all people.



HYHJ-1330 I

T1



HYHJ-1330 II

T2



HYHJ-1330 III

T3

Parameter

Monitor	3" color display
Camera	Pixel: ≥ 2 million; illuminance: ≥ 3000 lux
Battery	Battery life: ≥ 6 h; charging time: ≤ 4 h; capacity: 3400mAh
Function	Take pictures, video, freeze, store and export data
Storage	-40°C~70°C ambient temperature
Humidity	$\leq 93\%$ rh



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HYHJ-KM **Disposable** **Video Laryngoscope** /HYHJ-KM3C, HYHJ-KM5C/



Application

It is provided for clinical use to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and implement mechanical ventilation.

Suitable for anesthesia departments, ICU, pediatrics, emergency centers, ambulances,



Feature and Configuration

This product is divided into two models: HYHJ-KM3C and HYHJ-KM5C, has a total of 8 types of disposable laryngoscope blades (Miller00, Miller0, Miller1, Mac1, Mac2, Mac3, Mac4, Mac5), can be applied to all people.

Single-use laryngoscope blade



Parameter

Monitor	3" LCD display / 3.5" color touch screen display
Camera	pixel:≥3million; light illuminance:≥3000 lux
Battery	battery life:≥4h;charging time:≤5h;capacity:3400mAh
Interface Function	HDMI,Type-C;Photo,video,freeze,playback,store and export data
Storage	-40°C~70°C ambient temperature
Humidity	≤93% rh



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HYHJ-KP Disposable Video Laryngoscope /HYHJ-KP3C, HYHJ-KP5C/



Application

It is provided for clinical use to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and implement mechanical ventilation.

Suitable for anesthesia departments, ICU, pediatrics, emergency centers, ambulances,



Feature and Configuration

The monitor and camera parts of the product adopt plug-and-play connection mode, which is highly efficient, sturdy and beautiful; The product is divided into HYHJ-KP3C and HYHJ-KP5C two models, respectively corresponding to Miller00, Miller0, Miller1 and Mac1, Mac2, Mac3, Mac4, Mac5, there are 8 kinds of disposable laryngoscope blades, can be applied to all people.

Single-use laryngoscope blade



Parameter

Monitor	3" LCD display / 3.5" color touch screen display
Camera	pixel:≥3million; light illuminance:≥3000 lux
Battery	battery life:≥4h; charging time:≤5h; capacity:3400mAh
Interface Function	HDMI, Type-C; Photo, video, freeze, playback, store and export data
Storage	-40°C~70°C ambient temperature
Humidity	≤93% rh



Shenzhen Haiye Medical Technology Co., Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HY-YG5.2 Video Stylet



Application

It is provided for clinical use to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and implement mechanical ventilation.


Suitable for anesthesia departments, ICU, pediatrics, emergency centers, ambulances, etc. in medical institutions



Features and Configuration

Product structure: handle and insertion part

This product can be used with 3.5" touch screen display or 12.1" (non-touch), 13.3" big touch screen monitor.

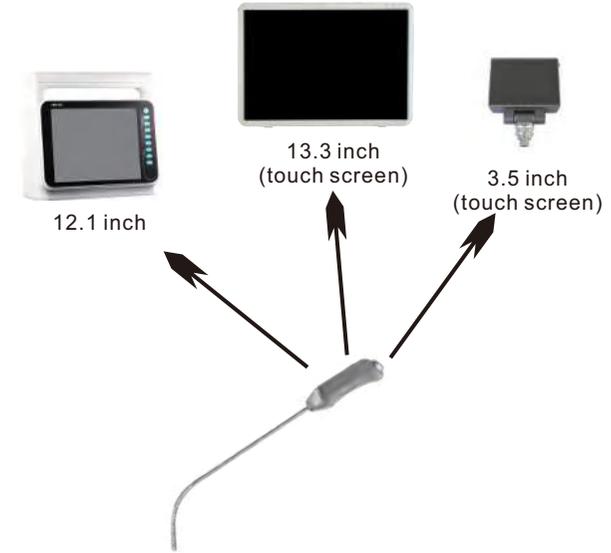

Field angle	70°±9°
Feild depth	3~100mm
Size	diameter 5.2±0.1mm, length 400±10mm, bending angle 110°±10°
Camera	2 million pixels; resolution ≥7.13lp/mm
Illumination	≥3000lx
Function	connect to 3.5", 12.1" or 13.3" display, one-click photo and video



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE




HP Series Reusable Anesthesia Fiber Optic Laryngoscope HP-W/HP-Z

Application

It is used clinically to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and mechanical ventilation.

Suitable for anesthesia departments, ICU, pediatrics, emergency centers, ambulances, etc. in medical institutions.

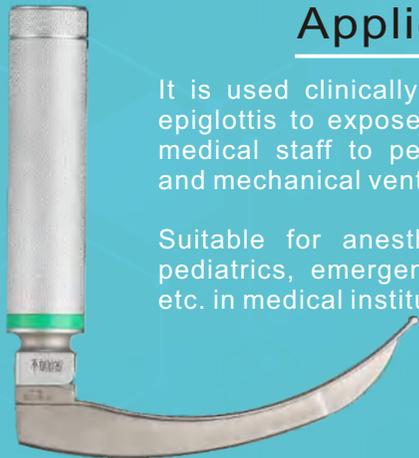

Feature

This product is divided into two series, HP-W and HP-Z, corresponding to adults and children respectively. Using high brightness LED cold light bulb, 4mm fiber bundle, high brightness, wide vision, built-in fiber optic design to avoid breakage and wear, convenient for cleaning and disinfection, and extend the service life; The laryngoscope blade is made of stainless steel and cast in one piece, which is made in one-step forming, using a three-point fixing method.



Parameter

Model of laryngoscope blade	Composition	Size of laryngoscope blade	Applicable user
HP-W1	Mac 4	155*15 mm	Adult
HP-W2	Mac 3	134*14.2 mm	Adult
HP-W3	Mac 2	114*14 mm	Adult
HP-Z1	Mac 1	93*10.9 mm	Pediatric
HP-Z2	Miller 1	80*9.4 mm	Pediatric
HP-Z3	Miller 0	68*10.7 mm	Pediatric




Shenzhen Haiye Medical Technology Co.,Ltd.
Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China
Tel: +86-29-87872952
Email: info@haiyemedical.com
Website: www.hyendoscope.com

HAIYE



HP-Y

Single-use Anesthesia Fiber Optic Laryngoscope

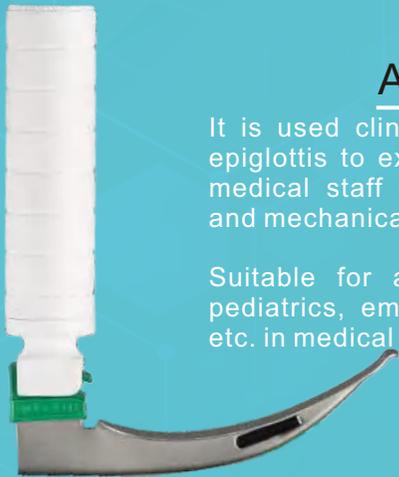

Application

It is used clinically to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and mechanical ventilation.

Suitable for anesthesia departments, ICU, pediatrics, emergency centers, ambulances, etc. in medical institutions.



Feature and Configuration

The entire product is packaged in sterile packaging and is for one-time use, requiring no cleaning or disinfection; The laryngoscope blade is made of stainless steel and is available in 6 specifications, which can meet the needs of different groups of people; The handle has a built-in high-brightness LED cold light source bulb, which is transmitted to the front end through the optical fiber bundle on the laryngoscope blade. High brightness, wide field of view.



Parameter

Model of laryngoscope blade	Shape of laryngoscope blade	Size of laryngoscope blade	Applicable user
Y1	Curved	153*13.8 mm	Adult
Y2	Curved	132*14.8 mm	Adult
Y3	Curved	100*13.8 mm	Adult
Y4	Straight	101.7*9.7 mm	Pediatric
Y5	Straight	80*10.2 mm	Pediatric
Y6	Straight	67*10.5 mm	Pediatric



Shenzhen Haiye Medical Technology Co., Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HYHJ-KX Series

Disposable Visual Laryngoscope Blade



Models:

HYHJ-KX00a, HYHJ-KX00b
HYHJ-KX0a, HYHJ-KX0b
HYHJ-KX1a, HYHJ-KX1b
HYHJ-KX2a, HYHJ-KX2b
HYHJ-KX3a, HYHJ-KX3b
HYHJ-KX4a, HYHJ-KX4b

Application

For clinical provocation of the patient's epiglottis to expose the vocal folds, guiding healthcare professionals to extubate and implement mechanical ventilation under visualization.

Applicable to the anesthesiology department, ICU, pediatrics department of medical institutions, emergency center, ambulance, etc.



Feature

The product is aseptically packaged for single use, minimizing the risk of cross infection.

The product is available in 6 sizes, and work with our video processors, it can be used by most people, from pediatric to adult.



Separate cable
and laryngoscope blade



one piece
(cable and laryngoscope blade)

Parameter

Monitor	12.1" or 13.3" touch screen (purchase separately)
Camera	one million pixels
Function	Take picture, record, freeze, playback, store and export data
Storage	ambient environment -40°C~70°C
Humidity	≤93% rh



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HY-ZQ series Electronic Bronchoscope

/HY-ZQ2.0, HY-ZQ2.8, HY-ZQ4.5, HY-ZQ5. 2, HY-ZQ5.8/



Feature and Configuration

Product production adheres to the spirit of craftsmanship, using sufficient materials and fine workmanship, and using a smaller outer diameter to achieve a larger instrument channel;

Product structure: 3 sizes monitors are optional, 3.5 inch, 12.1 inch and 13.3 inch are touchable screen, the tip part, the operating part and the insertion part; The products are divided into 4 specifications: HY-ZQ2.0, HY-ZQ2.8, HY-ZQ4.5, HY-ZQ5. 2 and HY-ZQ5.8.



3.5 inch
touch screen
image processor



13.3 inch
touch screen
image processor



HD3b
image processor

Parameter

Field angle	120°, tolerance ±10%
Field depth	3mm~150mm
Rotation	up 0° to 180°, down 0° to 130°, tolerance ±10%
Working time	≥240 mins
Instrument channel	1.2mm~3.0mm (≥Ø1.0mm, HY-ZQ2.0 not compatible)
Length of insert part	600mm, tolerance ±5%



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HY-BYH series Electronic Nasopharyngo Laryngoscope

HY-BYH1, HY-BYH2, HY-BYH3, HY-BYH4



Feature and Configuration

Product production adheres to the spirit of craftsmanship, using sufficient materials and fine workmanship, and using a smaller outer diameter to achieve a larger instrument channel;

Product structure: 2 sizes monitors are optional, 3.5 inch and 13.3 inch are touchable screen. HD3b model image processor can be connected to high-definition medical displays for better display effects.; The products are divided into 4 specifications: HY-BYH1, HY-BYH2, HY-BYH3, HY-BYH4



3.5 inch
touch screen
image processor



13.3 inch
touch screen
image processor



HD3b
image processor

Parameter

Field angle	120°, tolerance ±10%
Field depth	3mm~150mm
Rotation	up 0° to 130°, down 0° to 130°, tolerance ±10%
Working time	≥240 mins
Instrument channel	1.2mm~2.6mm (≥Ø1.0mm, HY-BYH1 not compatible)
Length of insert part	360mm, tolerance ±5%



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952


Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HY-MN2.5 Electronic Ureteroscope

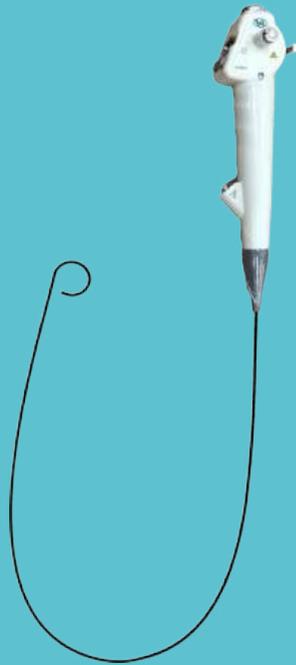

Feature and Configuration

Product production adheres to the spirit of craftsmanship, using sufficient materials and fine workmanship, and using a smaller outer diameter to achieve a larger instrument channel;

Product structure: 2 sizes monitors are optional, 3.5 inch and 13.3 inch are touchable screen. The HD3b model image processor can be connected to high-definition medical displays for better display effects.



3.5 inch
touch screen
image processor



13.3 inch
touch screen
image processor



HD3b
image processor

Parameter

Field angle	120°, tolerance $\pm 10\%$
Field depth	3mm~150mm
Rotation	up / down 275° , tolerance 10% \pm
OD of insert part	2.5mm
ID of insert part	1.2mm
Length of insert part	700mm, tolerance $\pm 5\%$



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China
Tel: +86-29-87872952
Email: info@haiyemedical.com
Website: www.hyendoscope.com

HAIYE



HY-YM2.8

Single-use Flexible Electronic Ureteroscope

HY-YM2.8A, HY-YM2.8B



Application

Used in conjunction with the electronic endoscope image processor produced by our company, it can be inserted into the urethra and renal pelvis through the urethra for observation, diagnosis, photography and treatment of the urethra and renal pelvis.



Configuration

The product is a disposable handheld flexible electronic endoscope consisting of an electronic ureteropyeloscope, a perfusion plug and an optional accessory waterproof cover.



Video converter



13.3 inch touch screen image processor



HD3b image processor

Parameter

Field angle	120°, tolerance ±15%
Field depth	3 ~100mm
Rotation	up/down 275°, tolerance ±15%
OD of insert part	2.5mm
ID of insert part	≥1.2mm
Length of insert part	700mm



Shenzhen Haiye Medical Technology Co., Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HY-YPG

Single-use Flexible Electronic Cystoscope

HY-YPG3.9, HY-YPG-4.5, HY-YPG5.3



Application

Used in conjunction with the electronic endoscope image processor produced by our company. It is inserted into the urethra and renal pelvis through the urethra and is used for observation, diagnosis, photography and treatment of the urethra and renal pelvis.



Configuration

This product is a disposable handheld flexible electronic endoscope. The operating part is designed in accordance with ergonomics. It can take pictures, record videos, and switch freely. The whole machine is a disposable design, which effectively prevents cross infection.



Video converter



13.3 inch touch screen image processor



HD3b image processor

Field angle	120°, tolerance ±10%
Field depth	3~100 mm
Camera pixels	400*400
Product specifications	HY-YPG3.9; HY-YPG4.5; HY-YPG5.3
ID of insert part	≥1.4-2.8 mm ± 0.2 mm
OD of insert part	≥3.9-5.3 mm ± 0.3 mm



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HY-YZQ series

Disposable Bronchoscope

/HY-YZQ2.0, HY-YZQ2.5, HY-YZQ4.5, HY-YZQ5.3/



Usage

Used in conjunction with the electronic endoscope image processor produced by our company, it is used by medical units to provide the necessary images for observation, diagnosis, photography and bronchi and treatment of the trachea and bronchi of adults and children over 3 months old. Product structure: The product is mainly composed of bending part, inserting part, operation part, conical sheath and video connection line.



Feature and Configuration

Product production adheres to the spirit of craftsmanship, using sufficient materials and fine workmanship, and using a smaller outer diameter to achieve a larger instrument channel;

Product structure: 3 sizes monitors are optional, decoder box, 12.1 inch and 13.3 inch are touchable screen, the tip part, the operating part and the insertion part; The products are divided into 4 specifications: HY-YZQ2.0, HY-YZQ2.5, HY-YZQ4.5 and HY-YZQ5.3.



Video converter



13.3 inch touch screen image processor



HD3b image processor

Parameter

Field angle	120°, tolerance ±10%
Field depth	3 ~100mm
Rotation	up 0° to 180°, down 0° to 130°
Camera	CMOS electronic imaging
Instrument channel	1.2mm~ 2.8mm (≥Ø1.0mm, HY-YZQ2.0 not compatible)
Length of insert part	600mm, tolerance ±5%



Shenzhen Haiye Medical Technology Co., Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HY-YBYH series Single-use Nasopharyngoscope

/HY-YBYH2.0, HY-YBYH2.5, HY-YBYH3.9, HY-YBYH4.5, HY-YBYH5.3 /



Feature and Configuration

Product production adheres to the spirit of craftsmanship, using sufficient materials and fine workmanship, and using a smaller outer diameter to achieve a larger instrument channel;

Product structure: 2 sizes monitors are optional, 3.5 inch and 13.3 inch are touchable screen. The HD3b model image processor can be connected to high-definition medical displays for better display effects.

The products are divided into 5 specifications:

HY-YBYH2.0, HY-YBYH2.5, HY-YBYH3.9, HY-YBYH4.5, HY-YBYH5.3.



Video converter



13.3 inch
touch screen
image processor



HD3b
image processor

Parameter

Field angle	120°, tolerance ±10%
Field depth	3mm~100mm
Rotation	up 0° to 210°, down 0° to 210, tolerance ±10%
Working time	≥240 mins
Instrument channel	1.2mm~2.8mm (≥Ø1.0mm, HY-YBYH2.0 not compatible)
Length of insert part	380mm ,tolerance ±5%



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HY-SP1 Video Converter




Application

It can be compatible with most of the endoscopes produced by Haiye, converting video footage to smartphones and computers in real time via OTG adapter cables and Type-C data transfer cables.



Feature

The Video Converter (decoder) supporting the conversion of video images of the following endoscopes of Haiye, the corresponding models and product names include: HYHJ-KX disposable video laryngoscope, HYHJ-KC reusable video laryngoscope blade, HYHJ-KP camera system and disposable blades, HY-YG5.2 video stylet, HY-ZQ reusable bronchoscope, HY-YMN disposable ureterscope, HY-YMN disposable bronchoscope, HY-YPG disposable cystoscope, as shown in the figure below:

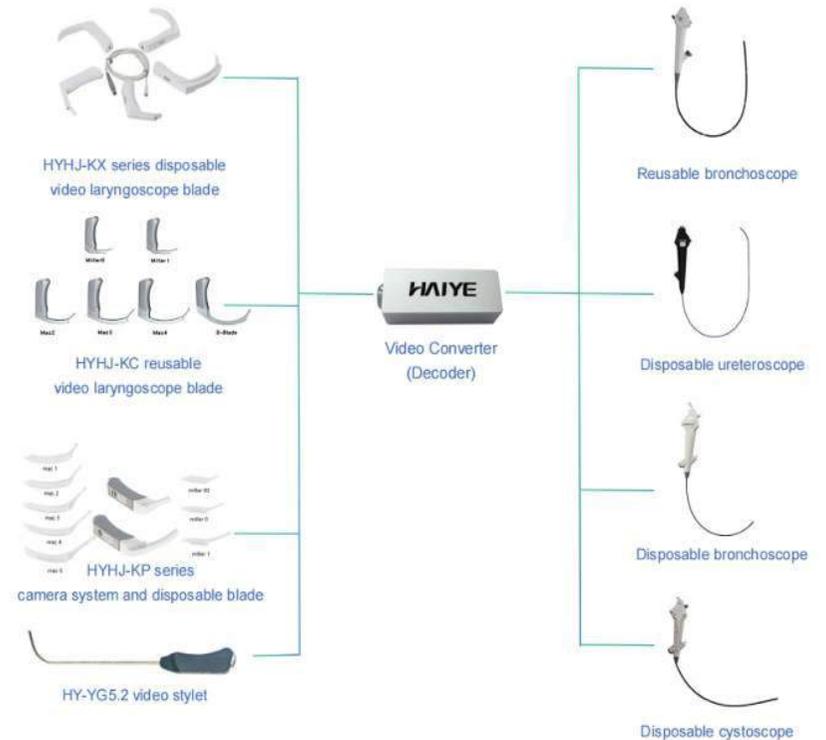

All disposable products: **No cross infection, No maintenance, No sterilization**



Shenzhen Haiye Medical Technology Co.,Ltd.
Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China
Tel: +86-29-87872952
Email: info@haiyemedical.com
Website: www.hyendoscope.com

HAIYE



HY-SCSTX1 Endoscopic Image Processor



Application

This product consists of an image processing section and is used in conjunction with our company's series of reusable endoscopes for internal surgical area image acquisition, and storage processing during endoscopic surgery.

This product is used in conjunction with the electronic endoscope produced by our company, to process the images captured by the electronic endoscope.

Feature

This product has complete functions, stable performance, and powerful image processing capabilities; Such as photography, video recording, storage, playback, white balance, etc. It is easy to carry and can seamlessly connect to our reusable endoscope.

Parameter

Built-in display	3.5", touch screen, 640*480
Compatible models	Reusable bronchoscope, Reusable nasopharyngoscope , Reusable ureteroscope
Video output	HDMI, Type-C
Function	Take pictures, record, freeze, play back, white balance, etc
Language	Optional (English, Spanish, French, Arabic, Russian, Chinese, German, Portuguese, Japanese, Thai, Vietnamese, Korean)
Aspect ratio	4:3
Working temperature	-10°C~40°C
Working power	Voltage: a.c.220V; Frequency: 50Hz; Input power: 30VA
Memory function	Standard with 32GB memory card, storage format:jpg,mov.



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE




HY-TS13 Endoscopic Image Processor



Application

This product consists of an image processing section and is used in conjunction with our company's series of endoscopes for internal surgical area image acquisition, image enlargement, and storage processing during endoscopic surgery.

This product is used in conjunction with the electronic endoscope produced by our company, to process the images captured by the electronic endoscope.

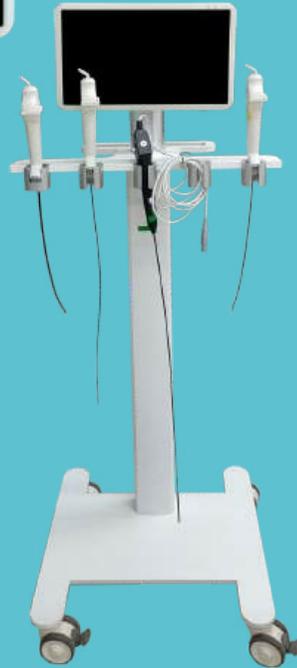

Feature

The product has complete functions, stable performance, and powerful image processing capabilities; such as photography, video recording, storage, playback, picture enlargement/reduction, white balance, etc. It has a built-in LED cold light source and multiple interfaces that can be seamlessly connected with our products.

Parameter

Built-in display	13.3", touch screen, 1920*1080
Compatible models	Video laryngoscope, bronchoscope, ureteroscope, cystoscope, nasopharyngoscope
Video output	DVI, HDMI, USB
Function	Take pictures, record, freeze, play back, zoom, white balance, etc
Language	Optional (English, Spanish, French, Arabic, Russian, Chinese, German, Portuguese, Japanese, Thai, Vietnamese, Korean)
Aspect ratio	4:3
Working temperature	5°C~40°C
Working power	Voltage: a.c.220V; Frequency: 50Hz; Input power: 30VA
Optional	Trolley, printer, computer, keyboard, etc



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HD3b Endoscope Image Processor



Application

It is suitable for medical institutions to connect with our endoscopes during clinical diagnosis or treatment to collect, process, data store and transfer the images of the human body cavity field of view observed by the endoscope, and transmit them to the display for output. The endoscope provides power for the user of the processor for related endoscopic surgeries; and can provide an illumination source for the endoscope.



Feature

The product has complete functions, stable performance, and powerful image processing capabilities; such as photography, video recording, storage, playback, picture enlargement/reduction, white balance, etc. It has a built-in LED cold light source and multiple interfaces that can be seamlessly connected with our products.

Parameter

Compatible models	Video laryngoscope, bronchoscope, ureteroscope, cystoscope, nasopharyngoscope
Video output	DVI, HDMI, USB
Function	Take pictures, record, freeze, play back, zoom, white balance, etc
Light source	Built-in cold light source
Size	374mm×340mm×123mm, tolerances±5%
Working temperature	5°C~40°C
Working power	Voltage: a.c.220V; Frequency: 50Hz; Input power: 30VA
Optional	Trolley, printer, computer, keyboard, etc



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HY-TC1 Mobile Stand




Application

Equipped with electronic endoscope image processor produced by our company, it can be used with our video laryngoscope, disposable uroscope, bronchoscope and other products for clinical tracheal intubation or treatment.

It is suitable for ENT, operating room, ICU, pediatrics, emergency center, ambulance, respiratory medicine, infection department, thoracic surgery, urology and so on.

Available Product

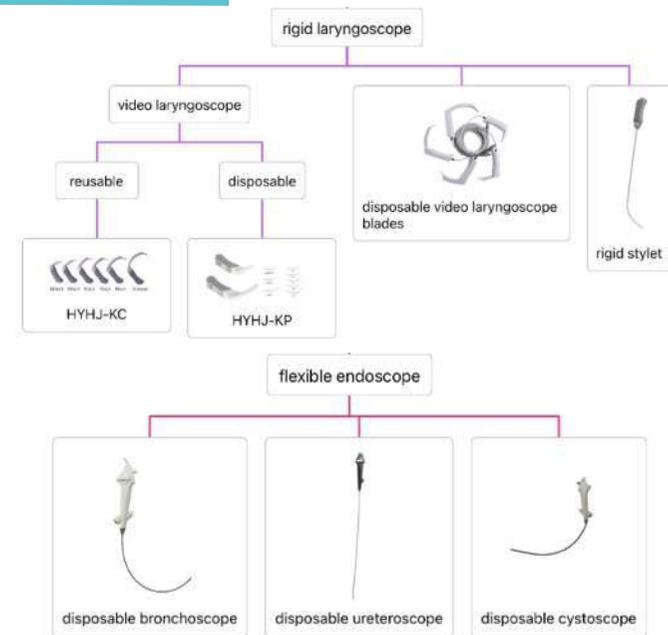

Product Description

The Mobile Stand consists of a support plate, support rods, workbench, multi-layer racks, bottom storage box, chassis, and universal casters. It can carry display screens of various specifications. The rack can hold image processors, endoscope camera systems, endoscope cold light sources, and printers, etc.



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952


Email: info@haiyemedical.com

Website: www.hyendoscope.com

HAIYE



HY-TC2 Multifuntional Airway Trolley

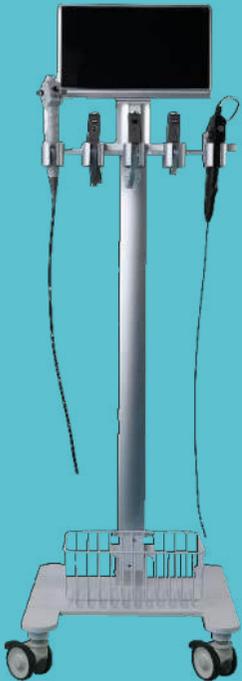

Application


Equipped with 12.1"HD1 electronic endoscope image processor produced by our company, it can be used with our video laryngoscope, disposable uroscope, bronchoscope and other products for clinical tracheal intubation or treatment.

It is suitable for ENT, operating room, ICU, pediatrics, emergency center, ambulance, respiratory medicine, infection department, thoracic surgery, urology and so on.

Parameter

Screen	12.1" / 13.3"color display
Resolution	1024*768 / 1920*1080
Field angle	89°
Function	HDMI, USB, signal , built-in storage card, support data storage and export, take photo, video,freeze frame
Compatible products	Video laryngoscope, video stylet, electronic bronchoscope, disposable ureteroscope, etc
Battery	6600mAh
Rated voltage	7.2V
Charging voltage	8.4V
Battery life	≥4h
Charging time	≤5h
Battery lifetime	500 cycles ≥80% of rated capacity
Accessories	adapter, power cord, data cable, manual

Scenario

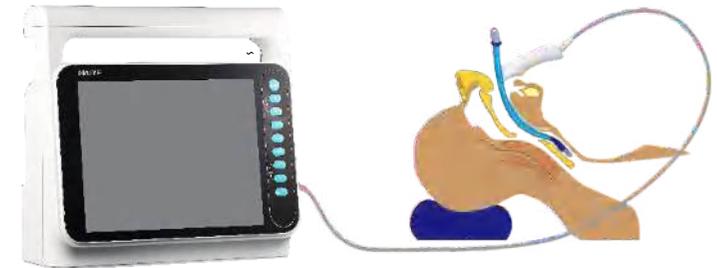

Shenzhen Haiye Medical Technology Co.,Ltd.


Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952


Email: info@haiyemedical.com

Website: www.hyendoscope.com

Haiye Full Range Video Laryngoscope

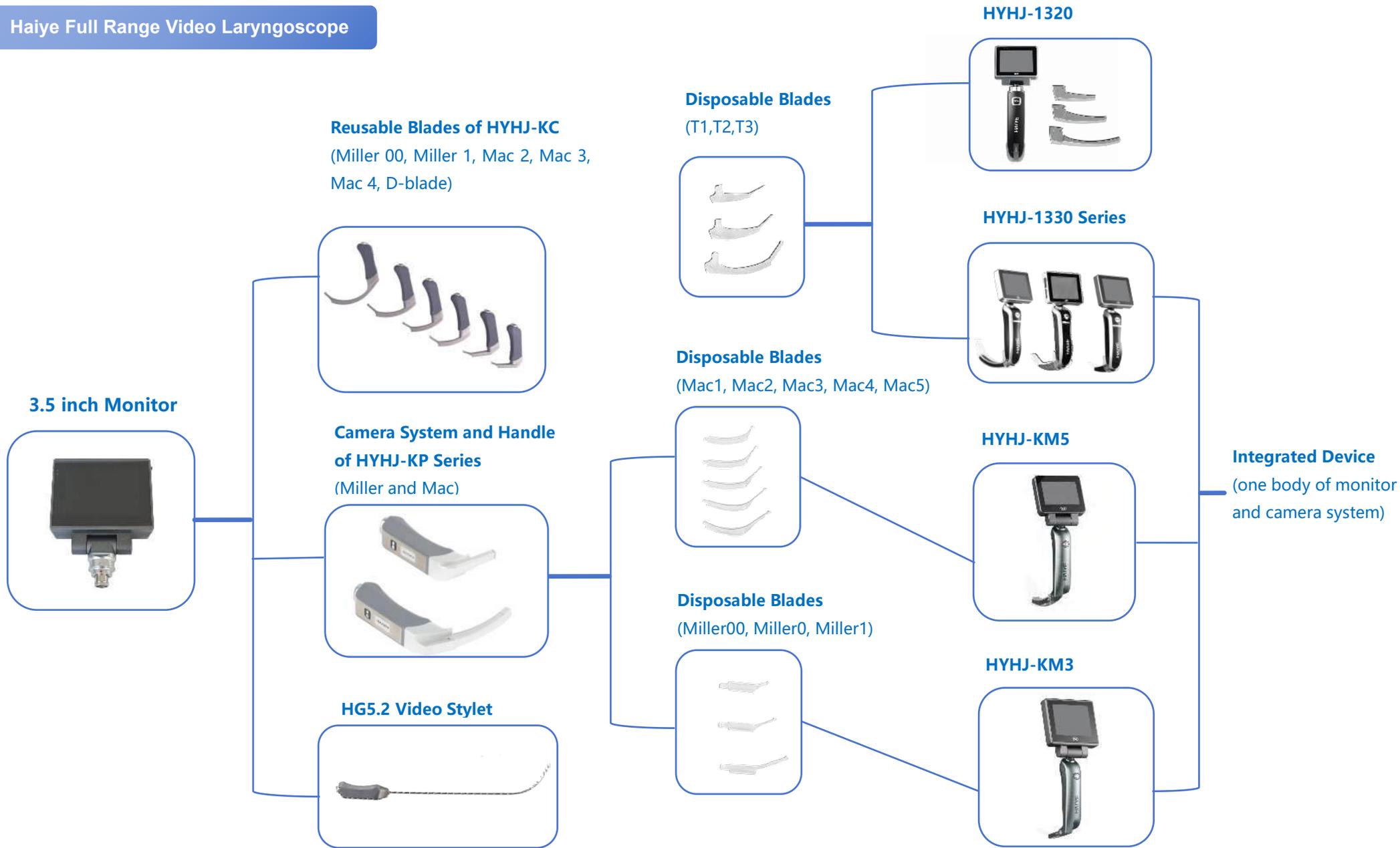

HAIYE Rigid Flexible Video Endoscope

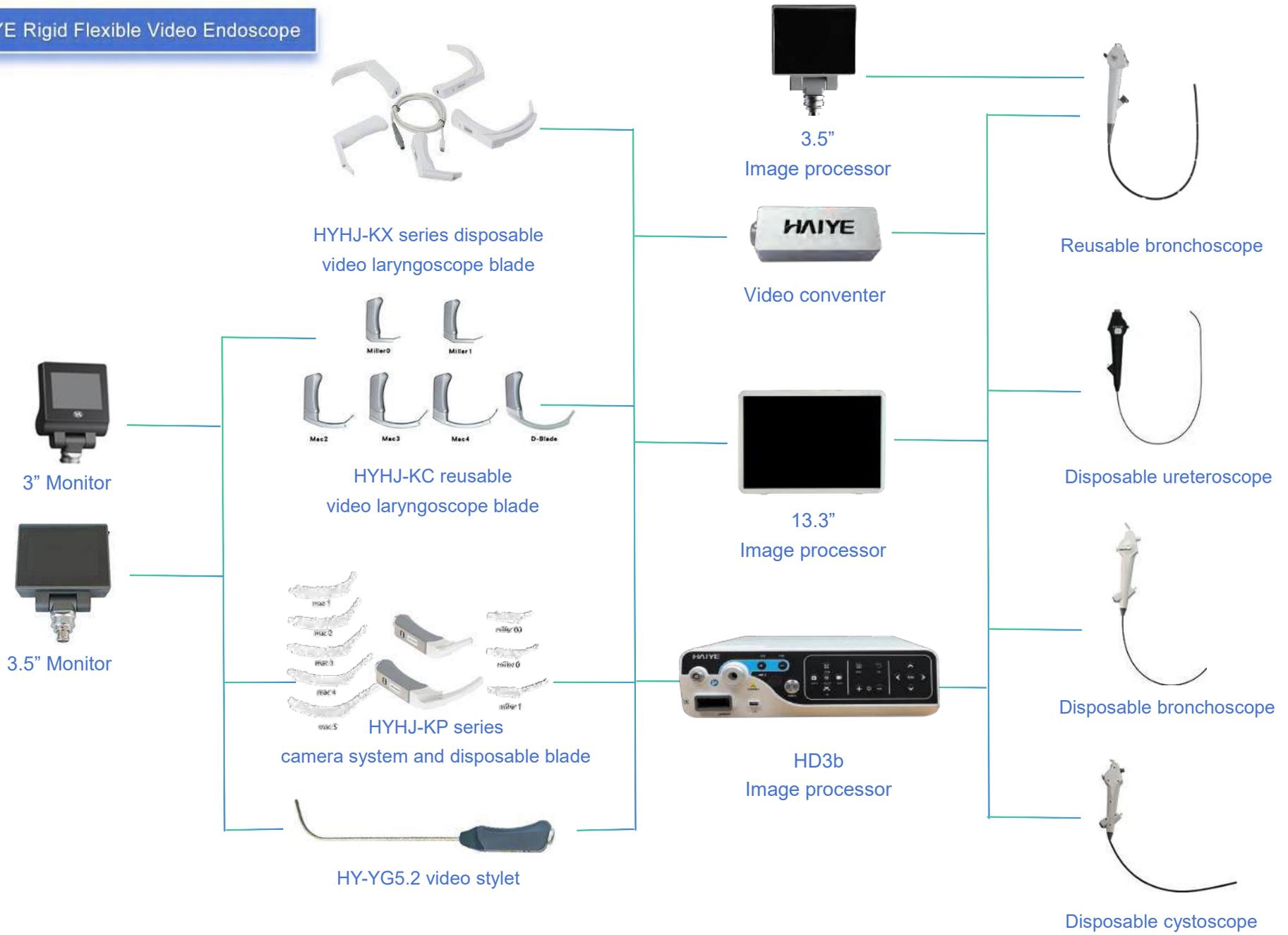

All disposable products: No cross infection, No maintenance, No sterilization