

HAIYE



HYHJ-KC Reusable Video Laryngoscope



Application


It is provided for clinical use to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and implement mechanical ventilation.

Suitable for anesthesia departments, ICU,



Feature and Configuration

The host display and camera parts are connected by plug and pull mode, which is efficient, strong and beautiful. The product can be used with 6 different specifications of blades can be reused after disinfection, its specifications are Miller 0, Miller 1, Mac 2, Mac 3, Mac 4, D-blade, can be applied to all people.

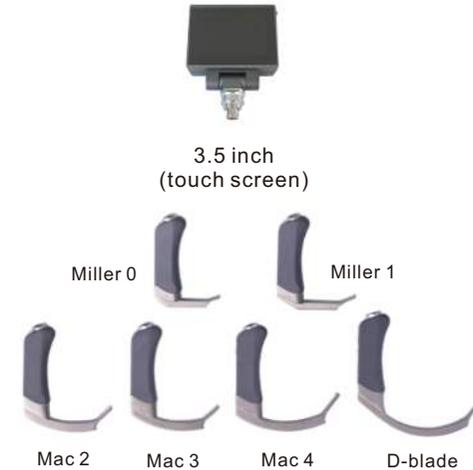

Parameter

Monitor	3" LCD display / 3.5" color touch screen display
Camera	pixel: ≥3million; light illuminance: ≥3000 lux
Battery	battery life: ≥4h; charging time: ≤5h; capacity: 3400mAh
Interface Function	HDMI, Type-C; Photo, video, freeze, playback, store and export data



Shenzhen Haiye Medical Technology Co.,Ltd.

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com